



PRODUCT DESCRIPTION

The LifeTec™ (P)-SRF V filter element is a pleated depth filter element with inner and outer liners and end caps made from stainless steel. The retention rate is >99.99999995 % related to 0.2 µm and 0.02 µm. The retention for nano-sized particles (2 nm – 10 nm) is larger than 99.99999991 % as verified in a DIN EN 1822 adopted test.

The LifeTec™ (P)-SRF V was developed for the safe and sterile filtration of air and other process gases in venting applications. The filter element fulfils the high requirements in food and beverage industries (dairies, food processing, soft drinks) as well as in pharmaceutical industries and works reliable even under extreme operating conditions.

The depth filter medium complies with the FDA (Food and Drug Administration) requirements CFR Title 21 (211.72 latest edition) and (EC) No. 1935/2004 for indirect food contact use.

The sturdy stainless-steel construction permits more than 160 possible sterilization cycles at specified conditions and withstands high differential pressures in both flow directions. LifeTec™ (P)-SRF V sterile filter elements ensure a safe and reproducible production.

MAIN FEATURES & BENEFITS

- Developed for the sterile filtration of air and gases in venting applications in which CIP reagents are applied
- High retention rate (bacteria, viruses and particles) down to 3 nm (nanometers) to ensure product and process integrity
- High temperature and mechanical resistance for outstanding performance, minimizes production down time and maintenance costs
- Suitable for sterilization, using VPHP and ozone - this leads to reduced total cost of ownership
- Resistant towards aerosolized CIP reagents

INDUSTRIES



- Dairies



- Food



- Breweries



- Pharmaceutical



- Chemical

PRODUCT SPECIFICATIONS

Features	Specifications
Retention rate	> 99.99999995 % at 0.2 µm > 99.99999995 % at 0.02 µm > 99.99999991 % at 0.003 µm
Filtration surface	0.3 m² per 10" element (254 mm) For other element sizes see correction factor CF in section "Available end cap configurations"
Operating Temperature	-20 °C (-4 °F) to +200 °C (+392 °F), at dry conditions
Maximum differential pressure	5 bar (73 psid) for -20 °C (-4 °F) up to +200 °C (392 °F), independent of the system pressure or flow direction
Typical compressed air service life time	12 months
Typical vent service life time	6 months
Sterilisation cycles*	121 °C (250 °F) saturated steam: >160 cycles (30 minutes) 131 °C (270 °F) saturated steam: >160 cycles (20 minutes) 141 °C (290 °F) saturated steam: >160 cycles (10 minutes) Independent of flow direction; forward and reversed steam flow possible
VPHP suitable	130 °C @ > 5000 ppm H ₂ O ₂ > 50 hours

* Figures are based on lab tests to evaluate steaming resistance. Filter elements need to be checked in actual use. Contact Donaldson for recommended autoclaving / steaming procedures.

APPLICATIONS

The pleated sterile depth filter LifeTec™ (P)-SRF V is designed and developed for the following applications:

Filtration of air and gases

- Tank Ventilation
- Container Aeration
- Carbon Dioxide
- Fermentation Air
- Technical Gases

QUALITY ASSURANCE

All products have been inspected and released by Quality Assurance as having met the following requirements:

- All filters are fabricated without the use of chemical additives.
- All sterile filters are integrity tested according to ASTM D 2986-91 and EN 1822 to verify compliance with established quality and design specifications and to assure consistent and reliable performance
- A factory certificate of compliance (2.2) according to DIN EN 10204 is available upon request.

MATERIAL DECLARATION (US & EU)

All components of the LifeTec™ (P)-SRF V filter cartridge are FDA listed for food contact use in the Code of Federal Regulations (CFR), Title 21. Donaldson Filtration Deutschland GmbH confirms that all materials used for the LifeTec™ (P)-SRF V elements meet regulatory and legislative requirements and guidelines for indirect food contact as detailed in European Regulation (EC) Number 1935/2004.

Filter materials		CFR Title 21
Filtermedium	Borosilicate	§ 177.2260
Impregnation	PTFE	§ 177.1550
Up- and downstream support	Stainless steel 1.4301	§ 211.65
Inner and outer liner	Stainless steel 1.4301	§ 211.65
Up- and downstream support media	PTFE	§ 177.1550
End caps	Stainless steel 1.4301	§ 211.65
Potting material	Silicone	§ 177.2600
O-rings	Silicone	§ 177.2600
Alternative o-rings	EPDM	§ 177.2600
	FEP over silicone	§ 177.1550
	FEP over viton	§ 177.1550

RETENTION OF MIKROORGANISMS

The LifeTec™ (P)-SRF V sterile filter elements were challenged with a specified bacteria, phage and particle aerosol for a defined time. Down stream analysis of the filtered air was done using impactor or impinger.

$$LRV = \log_{10} \left(\frac{\text{Number of organisms in challenge}}{\text{Number of organisms in filtrate}} \right)$$

- Brevundimonas Diminuta ($\geq 0.2 \mu\text{m}$):
LRV > 9
- MS2 Coliphagen ($\geq 0.02 \mu\text{m}$):
LRV > 9
- Nano-sized particles ($\geq 0.003 \mu\text{m}$):
LRV > 10

FLOW CHARACTERISTICS

Type LifeTec™ (P)-SRF V		Flow at 7 barg [m³/h]	
housing	element	nominal*	maximum
0006	03/10	60	90
0009	04/10	90	120
0012	04/20	120	180
0018	05/20	180	270
0027	05/25	270	360
0036	07/25	360	480
0048	07/30	480	720
0072	10/30	720	1080
0108	15/30	1080	1440
0144	20/30	1440	1920
0192	30/30	1920	2280
0288	30/50	2880	4320

*The given nominal flow rate in the table represents 100 % nominal air flow in the diagram.


Pressure [barg]	0	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
Correction Factor [-]	0.13	0.25	0.38	0.50	0.63	0.75	0.88	1.00	1.13	1.25	1.38	1.50	1.63	1.75	1.88	2.00	2.13

Nominal and maximum flow for other pressures can be calculated with the above correction factors.


$$\dot{V}_{corrected} = \frac{\dot{V}_{nominal}}{\text{Correction Factor}}$$


DE-WETTING CHARACTERISTICS

De-wetting characteristic of a SRF V 10/3 P7 after steaming at 1 bar (121°C) for 30 minutes. Flow is 140 Nm³/h at 2 bar absolute. Normal conditions are reached after ~ 20 seconds.

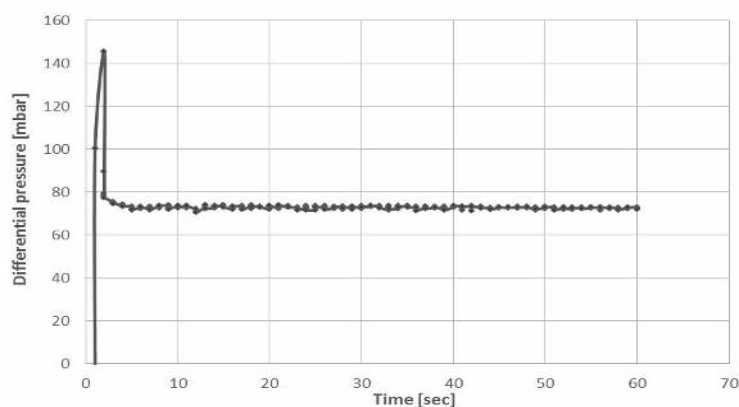

Nominal Flow P-SRF V 10/30 at 20°C

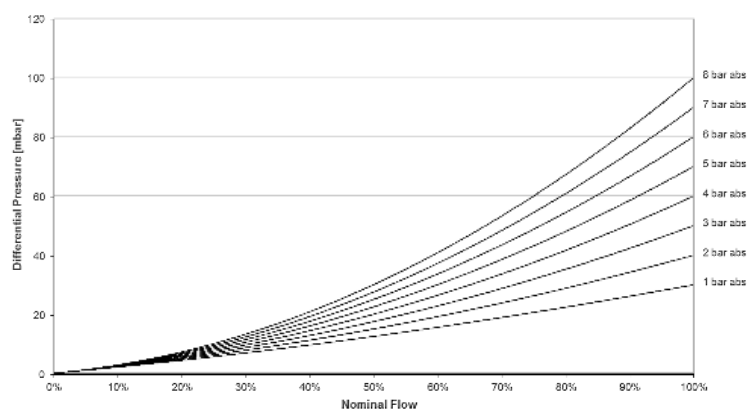

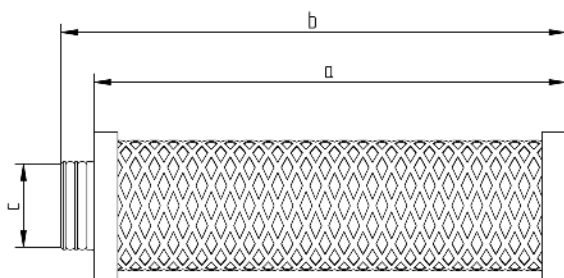


AVAILABLE END CAP CONFIGURATIONS

Dimensions of plug-connection							
Size	a		b		c*		CF**
	mm	inch	mm	inch	mm	inch	
03/10	76	2.99	87	3.42	30	1.18	0.15
04/10	104	4.09	118	4.64	30	1.18	0.20
04/20	104	4.09	118	4.64	37	1.46	0.20
05/20	128	5.04	142	5.59	37	1.46	0.25
05/25	128	5.04	142	5.59	37	1.46	0.34
07/25	180	7.08	194	7.64	37	1.46	0.49
05/30	128	5.04	142	5.59	61	2.40	0.49
07/30	180	7.08	196	7.71	61	2.40	0.70
10/30	254	10	270	10.63	61	2.40	1.00
15/30	381	15	397	15.63	61	2.40	1.51
20/30	510	20	526	20.63	61	2.40	2.02
30/30	764	30	780	30.63	61	2.40	3.03
30/50	764	30	780	30.63	89	3.50	3.03

* of plug-connection with double-o-ring

** Correction factor filtration surface

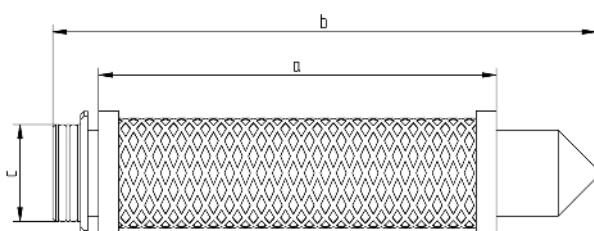

Dimensions Code 7 connection						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
5"	125	4.92	190	7.48	56,5	2.22
10"	250	9.84	315	12.40	56,5	2.22
20"	500	19.68	585	22.24	56,5	2.22
30"	750	29.53	815	32.08	56,5	2.22

CODE 7: 2 x 226 o-rings, bayonet 2 locking tabs, locating fin

OTHER END CAP CONFIGURATIONS ON REQUEST

For more information, please contact your Donaldson Sales Engineer and visit our website at www.donaldson.com!



donaldson.com/process

Donaldson Company, Inc.
Minneapolis, MN

Contact us



Important Notice: Many factors beyond the control of Donaldson can affect the use and performance of Donaldson products in a particular application, including the conditions under which the product is used. Since these factors are uniquely within the user's knowledge and control, it is essential the user evaluate the products to determine whether the product is fit for the particular purpose and suitable for the user's application. All products, specifications, availability and data are subject to change without notice, and may vary by region or country.

F119045 ENG (03/25) LifeTec™ (P)-SRF V ©2025 Donaldson Co., Inc. Donaldson and the color blue are marks of Donaldson Company, Inc. All other marks belong to their respective owners.